


iSLEEVE™
 Expandable Introducer Set

Included with the ACURATE neo™ Aortic Valve System
Ergonomic Hub

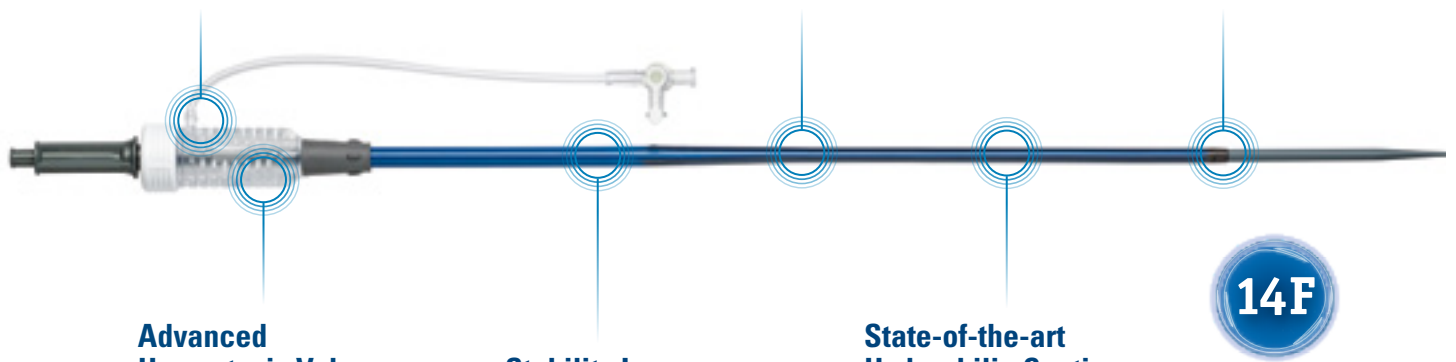
Provides total control

Responsive Tri-Fold Seams

Minimizes vessel trauma

Soft Radiopaque Tip

Enhances visualization


Advanced Hemostasis Valve

Minimizes blood loss

Stability Layer

Increases stability and sealing

State-of-the-art Hydrophilic Coating

Eases insertion and provides atraumatic removal

14F
iSLEEVE Introducer Set Specifications

Recommended Guidewire	0.035" (0.89 mm)
Sheath Working Length	32 cm
Dilator Length	48.7 cm
Valve Compatibility	For use with all ACURATE neo Aortic Valves
Use	Single Use
Sterilization	Ethylene Oxide Gas
Set Contents	<ul style="list-style-type: none"> • Introducer Sheath (1) • Dilator (1) • Product Mandrel (1) • Flushing Tube (1)
Carton Size	171.45 mm × 28.58 mm × 711.21 mm
Shelf Life	2 Years
Indications for Use	The 14F iSLEEVE Introducer Set is intended to facilitate femoral access to the vascular system for introduction and removal of the ACURATE™ TF Valve System and ancillary devices. The iSLEEVE Introducer Set is suitable for use in vessels ≥ 5.5 mm diameter.

Ordering Information

ACURATE <i>neo</i> [™] Aortic Valve System + 14F iSLEEVE [™] Expandable Introducer Set		
Ref/Catalog Number	Description	Units
BUN67ACURDSTFISL	ACURATE <i>neo</i> + 14F iSLEEVE Generic Bundle	1
BUN67ACURTFISL14FS*	ACURATE <i>neo</i> S + 14F iSLEEVE	1
BUN67ACUTFISL14FM*	ACURATE <i>neo</i> M + 14F iSLEEVE	1
BUN67ACUTFISL14FL*	ACURATE <i>neo</i> L + 14F iSLEEVE	1

* Only available in the following countries: Germany, Austria, Switzerland and the Netherlands.

All trademarks are the property of their respective owners. Illustrations for information purposes – not indicative of actual size or clinical outcome. All photographs taken by Boston Scientific.

CAUTION: The law restricts these devices to sale by or on the order of a physician. Indications, contraindications, warnings and instructions for use can be found in the product labeling supplied with each device. Information for use only in countries with applicable health authority registrations. Material not intended for use or distribution in France, Japan and the USA.

SH-526305-AE Printed in Germany by medicalvision.

**Boston
Scientific**
Advancing science for life[™]

www.bostonscientific.eu/ACURATEneo

© 2019 Boston Scientific Corporation
or its affiliates. All rights reserved.
DINSH0194EA